

# Cardiac Model Generator Software Detailed Design

Generated by Doxygen 1.12.0



<b>1 Namespace Index</b>	<b>1</b>
1.1 Package List	1
<b>2 Hierarchical Index</b>	<b>3</b>
2.1 Class Hierarchy	3
<b>3 Class Index</b>	<b>5</b>
3.1 Class List	5
<b>4 Namespace Documentation</b>	<b>7</b>
4.1 <code>_init_</code> Namespace Reference	7
4.1.1 Detailed Description	7
4.2 <code>CardiacModelGenerator</code> Namespace Reference	7
4.2.1 Detailed Description	8
4.2.2 Function Documentation	8
4.2.2.1 <code>clean_tetra_mesh()</code>	8
4.2.2.2 <code>generate_point_cloud()</code>	8
4.2.2.3 <code>generate_tetra_mesh()</code>	9
4.2.2.4 <code>get_cell_quality()</code>	9
4.2.2.5 <code>save_vista()</code>	10
<b>5 Class Documentation</b>	<b>11</b>
5.1 <code>CardiacModelGenerator.CardiacMeshalyzer</code> Class Reference	11
5.1.1 Detailed Description	13
5.1.2 Constructor & Destructor Documentation	13
5.1.2.1 <code>__init__()</code>	13
5.1.3 Member Function Documentation	14
5.1.3.1 <code>add_page()</code>	14
5.1.3.2 <code>clear_all_files()</code>	15
5.1.3.3 <code>close_program()</code>	15
5.1.3.4 <code>get_masks()</code>	16
5.1.3.5 <code>getMaskOverlay()</code>	16
5.1.3.6 <code>ImportDicomSeries()</code>	16
5.1.3.7 <code>on_clean_tetra_mesh()</code>	16
5.1.3.8 <code>on_extract_mesh_quality()</code>	17
5.1.3.9 <code>on_generate_point_cloud()</code>	18
5.1.3.10 <code>on_generate_tetra_mesh()</code>	18
5.1.3.11 <code>open_point_cloud_options()</code>	19
5.1.3.12 <code>save_point_cloud()</code>	19
5.1.3.13 <code>save_tetra_cloud()</code>	20
5.1.3.14 <code>setup_menu_bar()</code>	20
5.1.3.15 <code>show_page()</code>	21
5.2 <code>CardiacModelGenerator.CleanTetraMeshOptions</code> Class Reference	22
5.2.1 Detailed Description	23

5.2.2 Constructor & Destructor Documentation	24
5.2.2.1 <code>__init__()</code>	24
5.2.3 Member Function Documentation	25
5.2.3.1 <code>on_cancel()</code>	25
5.2.3.2 <code>on_ok()</code>	25
5.3 CardiacModelGenerator.HomePage Class Reference	26
5.3.1 Detailed Description	27
5.3.2 Constructor & Destructor Documentation	27
5.3.2.1 <code>__init__()</code>	27
5.3.3 Member Function Documentation	28
5.3.3.1 <code>load_dicom_series()</code>	28
5.3.3.2 <code>load_segmentation()</code>	29
5.3.3.3 <code>update_image()</code>	29
5.3.3.4 <code>view_set()</code>	30
5.4 CardiacModelGenerator.PointCloudOptions Class Reference	31
5.4.1 Detailed Description	32
5.4.2 Constructor & Destructor Documentation	32
5.4.2.1 <code>__init__()</code>	32
5.4.3 Member Function Documentation	33
5.4.3.1 <code>on_generate_point_cloud()</code>	33
5.4.3.2 <code>update_merging_tolerance_value()</code>	34
5.4.3.3 <code>update_point_size_value()</code>	34
5.4.4 Member Data Documentation	35
5.4.4.1 <code>merging_tolerance_slider</code>	35
5.5 CardiacModelGenerator.StartPage Class Reference	35
5.5.1 Detailed Description	36
5.5.2 Constructor & Destructor Documentation	37
5.5.2.1 <code>__init__()</code>	37
5.5.3 Member Function Documentation	38
5.5.3.1 <code>close_program()</code>	38
5.5.3.2 <code>load_image()</code>	38
5.5.3.3 <code>on_resize()</code>	39
5.5.3.4 <code>open_home_page()</code>	39
<b>Index</b>	<b>41</b>

# Chapter 1

## Namespace Index

### 1.1 Package List

Here are the packages with brief descriptions (if available):

<a href="#">_init_</a>	7
<a href="#">CardiacModelGenerator</a>	7



## Chapter 2

# Hierarchical Index

### 2.1 Class Hierarchy

This inheritance list is sorted roughly, but not completely, alphabetically:

wx.Dialog	
CardiacModelGenerator.CleanTetraMeshOptions . . . . .	<a href="#">22</a>
CardiacModelGenerator.PointCloudOptions . . . . .	<a href="#">31</a>
wx.Frame	
CardiacModelGenerator.CardiacMeshalyzer . . . . .	<a href="#">11</a>
wx.Panel	
CardiacModelGenerator.HomePage . . . . .	<a href="#">26</a>
wx.ScrolledWindow	
CardiacModelGenerator.StartPage . . . . .	<a href="#">35</a>





## Chapter 3

# Class Index

### 3.1 Class List

Here are the classes, structs, unions and interfaces with brief descriptions:

<a href="#">CardiacModelGenerator.CardiacMeshalyzer</a>	11
<a href="#">CardiacModelGenerator.CleanTetraMeshOptions</a>	22
<a href="#">CardiacModelGenerator.HomePage</a>	26
<a href="#">CardiacModelGenerator.PointCloudOptions</a>	31
<a href="#">CardiacModelGenerator.StartPage</a>	35



# Chapter 4

## Namespace Documentation

### 4.1 `_init_` Namespace Reference

#### 4.1.1 Detailed Description

Created on Thu Nov 7 20:05:18 2024

@author: vinayjani

### 4.2 CardiacModelGenerator Namespace Reference

#### Classes

- class [CardiacMeshalyzer](#)
- class [CleanTetraMeshOptions](#)
- class [HomePage](#)
- class [PointCloudOptions](#)
- class [StartPage](#)

#### Functions

- [generate\\_point\\_cloud](#) (coords1=None, masks1=None, coords2=None, masks2=None, coords3=None, masks3=None, whichmask=1, tol=0.1, colormap\_name="viridis", point\_size=5)
- [generate\\_tetra\\_mesh](#) (point\_cloud\_cleaned)
- [clean\\_tetra\\_mesh](#) (grid, subdivisions=2, poisson\_iterations=10, clean\_tolerance=0.001, quality↵ threshold=1e-5)
- [get\\_cell\\_quality](#) (final\_volumetric\_mesh)
- [save\\_vista](#) (filepath, pyvista\_object)
- **main** ()

#### Variables

- **BASE\_DIR** = os.path.abspath(os.path.join(os.path.dirname(\_\_file\_\_), ".."))
- str **IMAGE\_PATH** = BASE\_DIR+"/CardiacModelGenerator/static/mesh\_intro\_pic.png"

## 4.2.1 Detailed Description

Created on Tue Dec 3 17:35:28 2024

@author: vinayjani

## 4.2.2 Function Documentation

### 4.2.2.1 clean\_tetra\_mesh()

```
CardiacModelGenerator.clean_tetra_mesh (
    grid,
    subdivisions = 2,
    poisson_iterations = 10,
    clean_tolerance = 0.001,
    quality_threshold = 1e-5)
```

@brief Cleans and smooths a tetrahedral mesh to improve quality and smooth the mesh. This is done via a poisson subdivision of the tetrahedra comprising the mesh so as to improve resolution.

@param grid A PyVista UnstructuredGrid object representing the input volumetric mesh (a pv.core.pointset.UnstructuredGrid object).

@param subdivisions Number of subdivisions for mesh refinement (default: 2, int).

@param poisson\_iterations Number of smoothing iterations (default: 10, float).

@param clean\_tolerance Tolerance for cleaning the mesh (default: 0.001, float).

@param quality\_threshold Minimum acceptable quality for mesh cells (default: 1e-5, float).

@return A PyVista PolyData object representing the cleaned, smoothed, and increased resolution mesh (a pv.core.polydata.PolyData object).

Here is the caller graph for this function:



### 4.2.2.2 generate\_point\_cloud()


```
CardiacModelGenerator.generate_point_cloud (
    coords1 = None,
    masks1 = None,
    coords2 = None,
    masks2 = None,
    coords3 = None,
    masks3 = None,
    whichmask = 1,
    tol = 0.1,
    colormap_name = "viridis",
    point_size = 5)
```

```

@brief Generates and visualizes a point cloud from the segmentations and dicoms.
Atleast one coords and masks pair is needed.
@param coords1 Optional. First image stack coordinates (4D np.ndarray).
@param masks1 Optional. First set of masks corresponding to this view (a 3D np.ndarray).
@param coords2 Optional. Second image view coordinates(4D np.ndarray array).
@param masks2 Optional. Second set of masks corresponding to coords2 (a 3D np.ndarray) .
@param coords3 Optional. Third set of coordinates (4D np.ndarray array).
@param masks3 Optional. Third set of masks corresponding to coords3 (a 3D np.ndarray).
@param whichmask Mask value to extract points (default: 1, int).
@param tol Tolerance for cleaning the point cloud (default: 0.1, float).
@param colormap_name Colormap used for coloring the point cloud (default: "viridis", string which corresponds
@param point_size Size of the points in the visualization (default: 5, int).
@return A PyVista PolyData object representing the cleaned point cloud (a pv.core.pointset.PolyData object).

```

Here is the caller graph for this function:

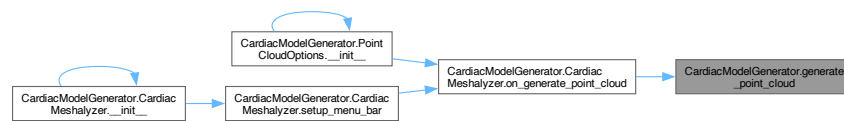

#### 4.2.2.3 generate\_tetra\_mesh()

```

CardiacModelGenerator.generate_tetra_mesh (
    point_cloud_cleaned)

```

```

@brief Generates a tetrahedral mesh from a cleaned point cloud. Note as long at the data is a pv.PolyData point
will work. It does not have to be cleaned in this tool the point cloud is cleaned.
@param point_cloud_cleaned A PyVista PolyData object representing the cleaned point cloud ( a pv.core.pointset
@return A PyVista UnstructuredGrid object representing the generated tetrahedral mesh (a pv.core.pointset.Unst

```

Here is the caller graph for this function:



#### 4.2.2.4 get\_cell\_quality()

```

CardiacModelGenerator.get_cell_quality (
    final_volumetric_mesh)

```

```

@brief Computes and visualizes cell quality for a tetrahedral mesh. This is done via computing the Jacobian of
of cell quality factors.
@param final_volumetric_mesh A PyVista UnstructuredGrid object representing the volumetric mesh (a pv.core.poi
@return A PyVista UnstructuredGrid object with cell quality values added as a scalar field (a pv.core.pointset

```

Here is the caller graph for this function:




#### 4.2.2.5 save\_vista()

```
CardiacModelGenerator.save_vista (  
    filepath,  
    pyvista_object)
```

@brief Saves a PyVista object to the specified file path.  
@param filepath The file path to save the PyVista object.  
@param pyvista\_object The PyVista object to be saved.

Here is the caller graph for this function:

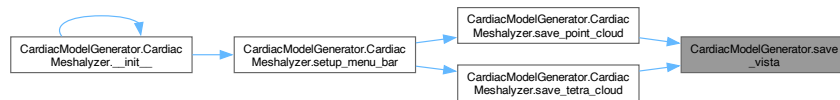


## Chapter 5


# Class Documentation

### 5.1 CardiacModelGenerator.CardiacMeshalyzer Class Reference

Inheritance diagram for CardiacModelGenerator.CardiacMeshalyzer:

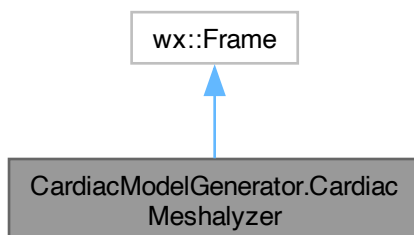

Collaboration diagram for CardiacModelGenerator.CardiacMeshalyzer:

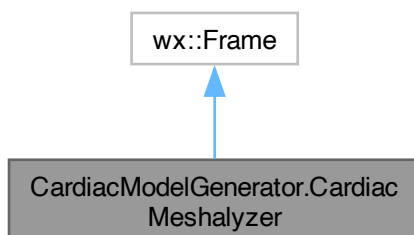

## Public Member Functions

- `__init__` (self, \*args, \*\*kwargs)
- `ImportDicomSeries` (self, folder\_path)
- `get_masks` (self, mask\_path)
- `setup_menu_bar` (self)
- `add_page` (self, name, page\_class)
- `show_page` (self, name)
- `save_point_cloud` (self, event)
- `save_tetra_cloud` (self, event)
- `clear_all_files` (self, event)
- `close_program` (self, event)
- `getMaskOverlay` (self, masks, volume)
- `open_point_cloud_options` (self, event)
- `on_generate_point_cloud` (self, event)
- `on_generate_tetra_mesh` (self, event)
- `on_clean_tetra_mesh` (self, event)
- `on_extract_mesh_quality` (self, event)

## Public Attributes

- dict `dicom_data` = {}
- dict `segmentation_data` = {}
- `panel` = wx.Panel(self)
- `sizer` = wx.BoxSizer(wx.VERTICAL)
- `page_container` = wx.Panel(self.panel)
- `page_sizer` = wx.BoxSizer(wx.VERTICAL)
- dict `pages` = {}
- dict `current_page` = None
- `on_generate_point_cloud`
- `on_generate_tetra_mesh`
- `on_clean_tetra_mesh`
- `on_extract_mesh_quality`
- `save_point_cloud`
- `save_tetra_cloud`
- `clear_all_files`
- `close_program`
- `last_point_cloud`
- `last_tetra_mesh`
- `colormap` = dialog.colormap
- `point_size` = dialog.point\_size
- `merging_tolerance` = dialog.merging\_tolerance
- `last_cleaned_mesh` = cleaned\_mesh
- `last_quality_mesh` = quality\_mesh



## 5.1.1 Detailed Description

```
@class CardiacMeshalyzer
@brief GUI application for managing and processing Dicoms and segmentations in the form of niftis. The
tool is most helpful for creating 3D tetrahedral meshes from this data
@details This class provides a graphical user interface (GUI) for handling DICOM series, generating point cloud
        creating tetrahedral meshes, and performing mesh cleaning and quality assessment.

@uml
@startuml
class CardiacMeshalyzer {
    - dicom_data : dict
    - segmentation_data : dict
    - panel : wx.Panel
    - sizer : wx.BoxSizer
    - page_container : wx.Panel
    - page_sizer : wx.BoxSizer
    - pages : dict
    - current_page : object

    + __init__(*args, **kwargs)
    + ImportDicomSeries(folder_path : str) : tuple
    + get_masks(mask_path : str) : numpy.ndarray
    + setup_menu_bar()
    + add_page(name : str, page_class : type)
    + show_page(name : str)
    + save_point_cloud(event : wx.Event)
    + save_tetra_cloud(event : wx.Event)
    + clear_all_files(event : wx.Event)
    + close_program(event : wx.Event)
    + getMaskOverlay(masks : numpy.ndarray, volume : numpy.ndarray) : numpy.ndarray
    + open_point_cloud_options(event : wx.Event)
    + on_generate_point_cloud(event : wx.Event)
    + on_generate_tetra_mesh(event : wx.Event)
    + on_clean_tetra_mesh(event : wx.Event)
    + on_extract_mesh_quality(event : wx.Event)
}
@enduml
```


## 5.1.2 Constructor & Destructor Documentation

### 5.1.2.1 \_\_init\_\_()


```
CardiacModelGenerator.CardiacMeshalyzer.__init__ (
    self,
    * args,
    ** kwargs)

@brief Initializes the CardiacMeshalyzer GUI application.
@param *args Positional arguments for the wx.Frame superclass.
@param **kwargs Keyword arguments for the wx.Frame superclass.
@details Sets up the GUI components, initializes global variables, and adds the Start Page and Home Page to th
```

Here is the call graph for this function:

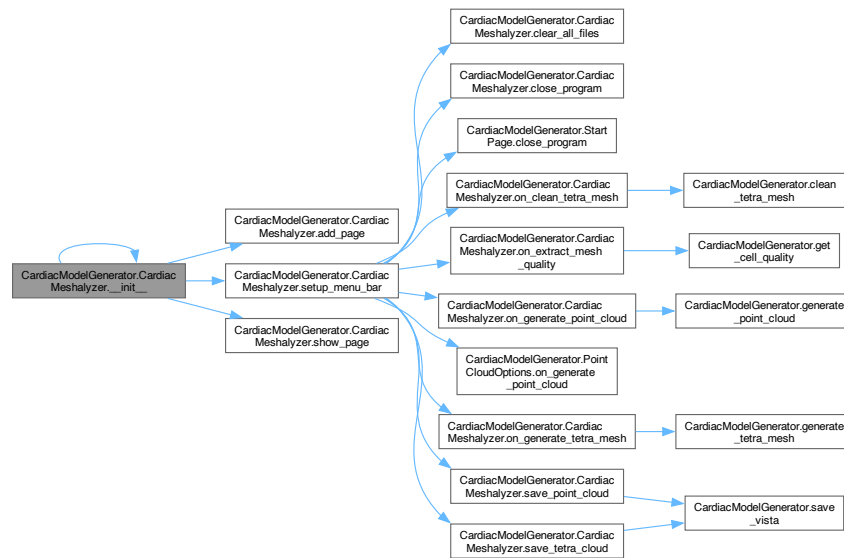

Here is the caller graph for this function:

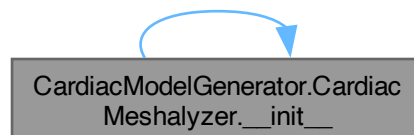

## 5.1.3 Member Function Documentation

### 5.1.3.1 add\_page()

```

CardiacModelGenerator.CardiacMeshalyzer.add_page (
    self,
    name,
    page_class)

```

@brief Adds a new page to the GUI.

@param name The name of the page to be added.

@param page\_class The class representing the page to add.

@details Creates an instance of the specified page class, adds it to the page container, and hides it initially.

Here is the caller graph for this function:




### 5.1.3.2 clear\_all\_files()

```
CardiacModelGenerator.CardiacMeshalyzer.clear_all_files (
    self,
    event)
```

`@brief` Clears all loaded DICOM data and segmentation files.  
`@param event` The wxPython event triggering this action.  
`@details` Resets the dictionaries holding DICOM and segmentation data and displays a confirmation message.

Here is the caller graph for this function:

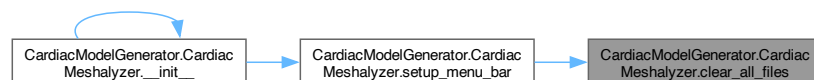


### 5.1.3.3 close\_program()

```
CardiacModelGenerator.CardiacMeshalyzer.close_program (
    self,
    event)
```

`@brief` Closes the application.  
`@param event` The wxPython event triggering this action.

Here is the caller graph for this function:

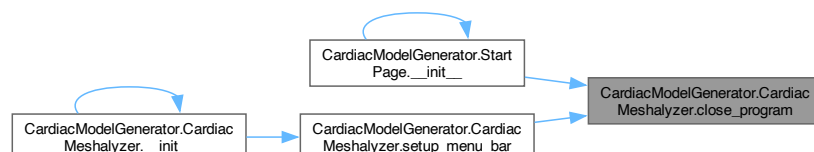

#### 5.1.3.4 get\_masks()

```
CardiacModelGenerator.CardiacMeshalyzer.get_masks (
    self,
    mask_path)
```

@brief Loads segmentation masks from a specified file.

@param mask\_path The file path to the segmentation mask in NIfTI format (string).

@return A NumPy array containing the segmentation mask data.

@details Reads the segmentation mask from the provided NIfTI file and converts it to a NumPy array for further

#### 5.1.3.5 getMaskOverlay()

```
CardiacModelGenerator.CardiacMeshalyzer.getMaskOverlay (
    self,
    masks,
    volume)
```

@brief Generates an overlay of the segmentation mask on the volume. The segmentation has some transparency.

@param masks A NumPy array containing the segmentation masks (a 3D np.ndarray).

@param volume A NumPy array representing the image volume (a 3D np.ndarray).

@return A NumPy array containing the overlay, where the segmentation masks are blended or burned with the volume.

@details This method normalizes the volume data, assigns colors based on the mask, and overlays colors on the volume to create a visual overlay. Handles up to 4 predefined segmentation classes and generates random colors.

@throws ValueError If the volume data is not numeric or not a NumPy array.

#### 5.1.3.6 ImportDicomSeries()

```
CardiacModelGenerator.CardiacMeshalyzer.ImportDicomSeries (
    self,
    folder_path)
```

@brief Processes a series of DICOM files in a given folder.

@param folder\_path The path to the folder containing DICOM files (string).

@return A tuple containing:

- Volume: A 3D NumPy array representing the reconstructed image volume.
- Coords: A 4D NumPy array with the absolute (x, y, z) coordinates for each pixel.
- out\_images: A NumPy array containing the original DICOM objects.

@details Reads DICOM files from the specified folder, reconstructs the image volume, and computes the absolute coordinates for each pixel using the image's orientation, spacing, and position metadata.

#### 5.1.3.7 on\_clean\_tetra\_mesh()


```
CardiacModelGenerator.CardiacMeshalyzer.on_clean_tetra_mesh (
    self,
    event)
```

@brief Handles the "Clean Tetra Mesh" menu option.


@param event The wxPython event triggering this action.

@details Opens a dialog to configure cleaning options and applies these settings to clean the tetrahedral mesh. Visualizes the cleaned mesh upon successful completion. Displays an error message if no tetrahedral mesh is found.

Here is the call graph for this function:

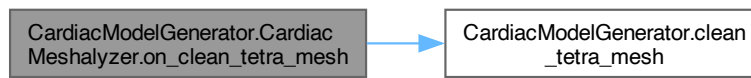

Here is the caller graph for this function:

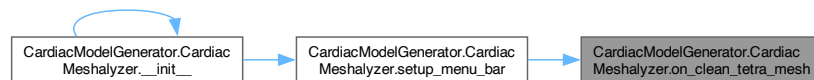

### 5.1.3.8 on\_extract\_mesh\_quality()


```
CardiacModelGenerator.CardiacMeshalyzer.on_extract_mesh_quality (
    self,
    event)
```

@brief Handles the "Extract Mesh Quality" menu option.


@param event The wxPython event triggering this action.

@details Computes the quality of cells in the tetrahedral mesh. Stores the quality mesh for future use and displays a success message upon completion. Displays an error message if no cleaned mesh exists.

Here is the call graph for this function:

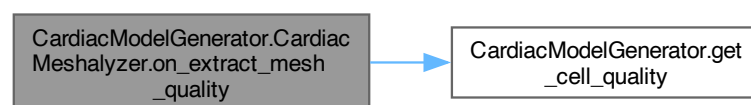

Here is the caller graph for this function:

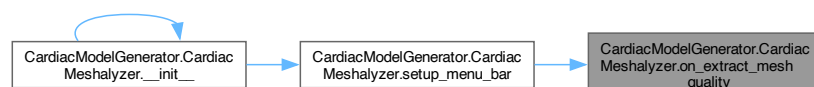

### 5.1.3.9 on\_generate\_point\_cloud()

```
CardiacModelGenerator.CardiacMeshalyzer.on_generate_point_cloud (
    self,
    event)
```

@brief Handles the "Generate Point Cloud" menu option.


@param event The wxPython event triggering this action.

@details Collects DICOM coordinate and mask data, opens the Point Cloud Options dialog for user input, and generates a point cloud using the specified options. If an error occurs during point cloud generation an error message is displayed.

Here is the call graph for this function:



Here is the caller graph for this function:

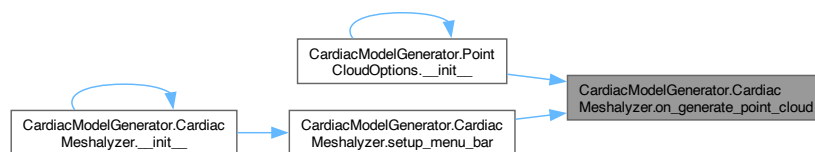

### 5.1.3.10 on\_generate\_tetra\_mesh()


```
CardiacModelGenerator.CardiacMeshalyzer.on_generate_tetra_mesh (
    self,
    event)
```

@brief Handles the "Generate Tetra Mesh" menu option.


@param event The wxPython event triggering this action.

@details Generates a tetrahedral mesh from the last generated point cloud. Displays an error message if no point cloud exists.

Here is the call graph for this function:

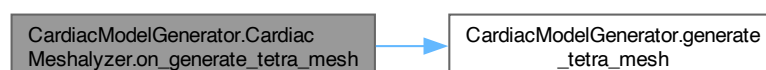

Here is the caller graph for this function:

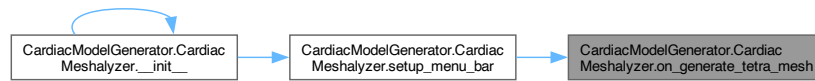

### 5.1.3.11 open\_point\_cloud\_options()

```
CardiacModelGenerator.CardiacMeshalyzer.open_point_cloud_options (
    self,
    event)
```

@brief Opens the Point Cloud Options dialog.

@param event The wxPython event triggering this action.

@details Displays a dialog to configure point cloud generation options, including colormap, point size, and mesh type. If the user confirms their selection, the 'generate\_point\_cloud' method is called with the selected options.

### 5.1.3.12 save\_point\_cloud()

```
CardiacModelGenerator.CardiacMeshalyzer.save_point_cloud (
    self,
    event)
```


@brief Saves the generated point cloud to a file using the save\_vista helper function.

@param event The wxPython event triggering this action.

Here is the call graph for this function:



Here is the caller graph for this function:

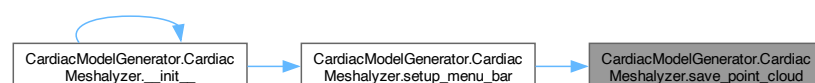

### 5.1.3.13 save\_tetra\_cloud()


```
CardiacModelGenerator.CardiacMeshalyzer.save_tetra_cloud (
    self,
    event)
```

@brief Saves the generated tetrahedral mesh to a file using the save\_vista helper function.  
 @param event The wxPython event triggering this action.

Here is the call graph for this function:



Here is the caller graph for this function:

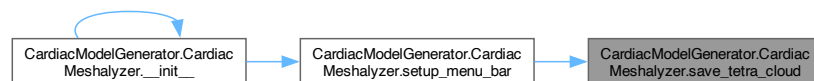


### 5.1.3.14 setup\_menu\_bar()

```
CardiacModelGenerator.CardiacMeshalyzer.setup_menu_bar (
    self)
```


@brief Configures and initializes the menu bar for the application.  
 @details Creates a menu bar with options for loading images, processing models, and saving/clearing data. Each menu item is bound to the corresponding event handler method.



Here is the call graph for this function:

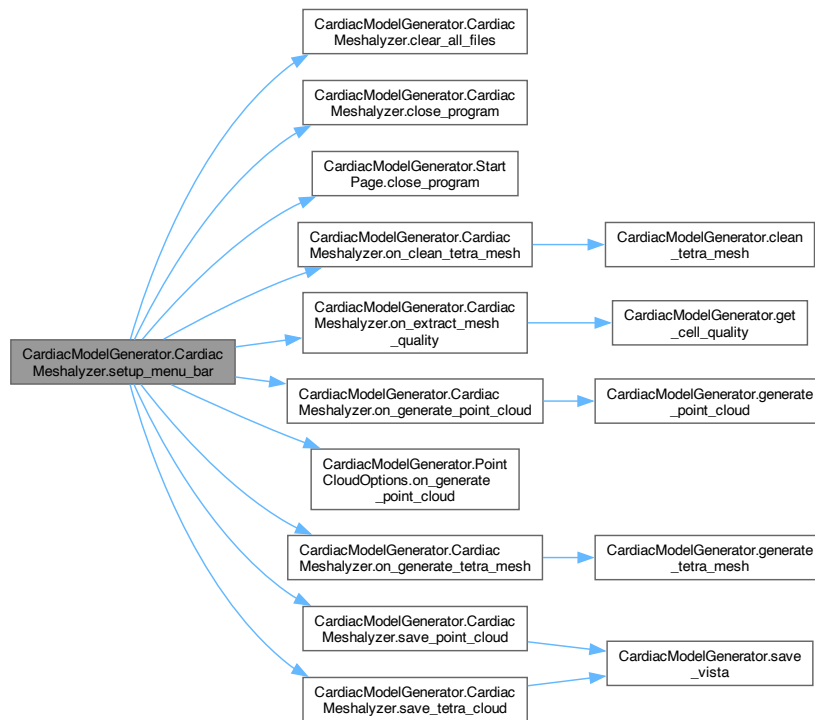

Here is the caller graph for this function:

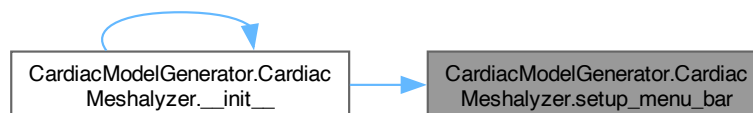

### 5.1.3.15 show\_page()

```
CardiacModelGenerator.CardiacMeshalyzer.show_page (
    self,
    name)
```

@brief Displays the specified page in the GUI.

@param name The name of the page to display.

@details Hides the currently visible page (if any) and shows the specified page.

Here is the caller graph for this function:




The documentation for this class was generated from the following file:


- `CardiacModelGenerator.py`

## 5.2 CardiacModelGenerator.CleanTetraMeshOptions Class Reference

Inheritance diagram for `CardiacModelGenerator.CleanTetraMeshOptions`:

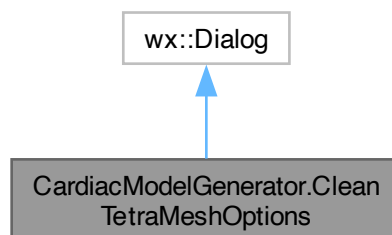

Collaboration diagram for `CardiacModelGenerator.CleanTetraMeshOptions`:

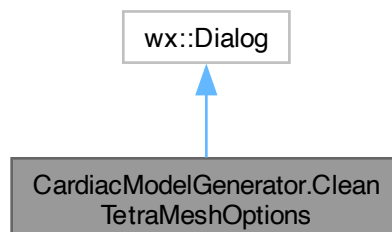

## Public Member Functions

- `__init__` (self, parent, \*args, \*\*kwargs)
- `on_ok` (self, event)
- `on_cancel` (self, event)

## Public Attributes

- `subdivisions_text` = wx.TextCtrl(self, value="2")
- `poisson_iterations_text` = wx.TextCtrl(self, value="10")
- `clean_tolerance_text` = wx.TextCtrl(self, value="0.001")
- `quality_threshold_text` = wx.TextCtrl(self, value="1e-5")
- `on_ok`
- `on_cancel`
- int `subdivisions` = 2
- int `poisson_iterations` = 10
- float `clean_tolerance` = 0.001
- int `quality_threshold` = 1e-5

### 5.2.1 Detailed Description

```
@class CleanTetraMeshOptions
@brief Dialog for configuring tetrahedral mesh cleaning options.
@details Provides controls for setting parameters such as subdivisions, Poisson iterations, cleaning tolerance
        and quality threshold for cleaning a tetrahedral mesh.

@uml
@startuml
class CleanTetraMeshOptions {
    - subdivisions_text : wx.TextCtrl
    - poisson_iterations_text : wx.TextCtrl
    - clean_tolerance_text : wx.TextCtrl
    - quality_threshold_text : wx.TextCtrl
    - subdivisions : int
    - poisson_iterations : int
    - clean_tolerance : float
    - quality_threshold : float

    + __init__(parent : wx.Window, *args, **kwargs)
    + on_ok(event : wx.Event)
    + on_cancel(event : wx.Event)
}

CleanTetraMeshOptions -- wx.Dialog : inherits
CleanTetraMeshOptions o-- wx.TextCtrl : "User input fields"
CleanTetraMeshOptions o-- wx.Button : "OK and Cancel buttons"
CleanTetraMeshOptions --> MeshCleaning : "Configures parameters for mesh cleaning"

' Notes for context
note top of CleanTetraMeshOptions
    CleanTetraMeshOptions allows users to define parameters for
    cleaning a tetrahedral mesh. It ensures proper numerical inputs
    and saves these configurations for further processing.
end note
@enduml
```


## 5.2.2 Constructor & Destructor Documentation

### 5.2.2.1 \_\_init\_\_()


```
CardiacModelGenerator.CleanTetraMeshOptions.__init__ (
    self,
    parent,
    * args,
    ** kwargs)
```

@brief Initializes the CleanTetraMeshOptions dialog.  
 @param parent The parent window that contains this dialog.  
 @param \*args Additional positional arguments for wx.Dialog.  
 @param \*\*kwargs Additional keyword arguments for wx.Dialog.  
 @details Creates a dialog with input fields for setting parameters related to tetrahedral mesh cleaning.  
 These parameters include:  
 - Subdivisions.  
 - Poisson Iterations.  
 - Cleaning Tolerance.  
 - Quality Threshold.  
 The dialog also provides OK and Cancel buttons for user interaction.

Here is the call graph for this function:

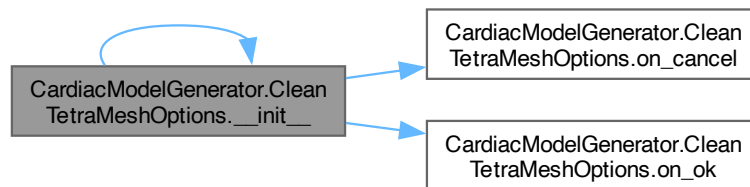

Here is the caller graph for this function:

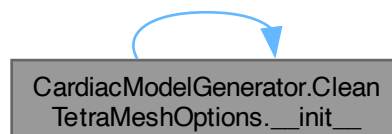

## 5.2.3 Member Function Documentation

### 5.2.3.1 on\_cancel()

```
CardiacModelGenerator.CleanTetraMeshOptions.on_cancel (  
    self,  
    event)
```

@brief Closes the dialog without saving user inputs.

@param event The wxPython event triggering this action.

@details Simply closes the dialog and returns a cancellation state without processing or saving any user input.

Here is the caller graph for this function:



### 5.2.3.2 on\_ok()

```
CardiacModelGenerator.CleanTetraMeshOptions.on_ok (  
    self,  
    event)
```


@brief Saves the user inputs and closes the dialog.

@param event The wxPython event triggering this action.

@details Reads and validates user inputs from the text fields for subdivisions, Poisson iterations, cleaning tolerance, and quality threshold. If the inputs are valid, they are saved as instance variables and the dialog is closed with a success state.

@throws ValueError If any input is not a valid numerical value.

Here is the caller graph for this function:

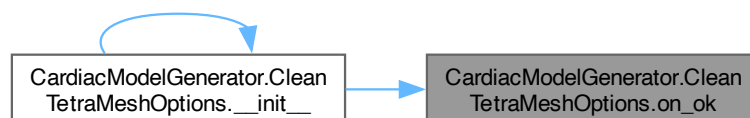


The documentation for this class was generated from the following file:


- CardiacModelGenerator.py

## 5.3 CardiacModelGenerator.HomePage Class Reference

Inheritance diagram for CardiacModelGenerator.HomePage:

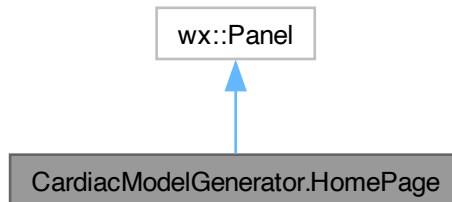

Collaboration diagram for CardiacModelGenerator.HomePage:

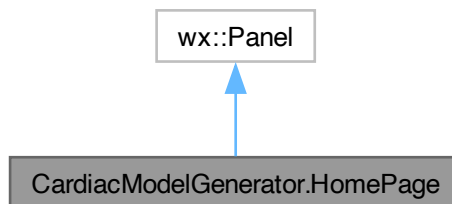

### Public Member Functions

- `__init__` (self, parent)
- `load_dicom_series` (self, view\_num)
- `load_segmentation` (self, view\_num)
- `view_set` (self, view\_num)
- `update_image` (self, event)

### Public Attributes

- `image_display` = `wx.StaticBitmap(right_panel)`
- `slider` = `wx.Slider(right_panel, minValue=0, maxValue=1, value=0, style=wx.SL_HORIZONTAL)`
- `update_image`
- `current_image_stack` = `None`

### 5.3.1 Detailed Description

```
@class HomePage
@brief Main page for interacting with DICOM images and segmentations.
@details This class provides functionality to load and display DICOM images and segmentations,
        manage image viewing with sliders, and bind buttons for different views and actions.

@uml
@startuml
class HomePage {
    - current_image_stack : np.ndarray

    + __init__(parent : wx.Window)
    + load_dicom_series(view_num : int)
    + load_segmentation(view_num : int)
    + view_set(view_num : int)
    + update_image(event : wx.Event)
}

HomePage *-- wx.Panel : inherits
HomePage o-- wx.Button : "Handles action buttons"
HomePage o-- wx.Slider : "Controls image stack navigation"
HomePage o-- wx.StaticBitmap : "Displays images"
HomePage --> DICOM : "Interacts with DICOM data"
HomePage --> Segmentation : "Handles segmentation data"

' Notes for additional context
note top of HomePage
    The HomePage class allows users to load, view, and interact
    with DICOM images and corresponding segmentation masks.
    It provides tools for navigating through image stacks and
    visualizing segmented regions overlaid on images.
end note
@enduml
```


### 5.3.2 Constructor & Destructor Documentation

#### 5.3.2.1 \_\_init\_\_()


```
CardiacModelGenerator.HomePage.__init__ (
    self,
    parent)

@brief Initializes the HomePage class and its UI components.
@param parent The parent wx.Window to which this HomePage belongs.
@details This method sets up the main layout of the HomePage, including:
    - A dark-themed background.
    - Buttons for loading image views and segmentations.
    - Buttons for viewing specific segmentations.
    - An image display area and a slider for navigating through image stacks.
    The layout is divided into a left panel for action buttons and a right panel for viewing and interact
```

Here is the call graph for this function:

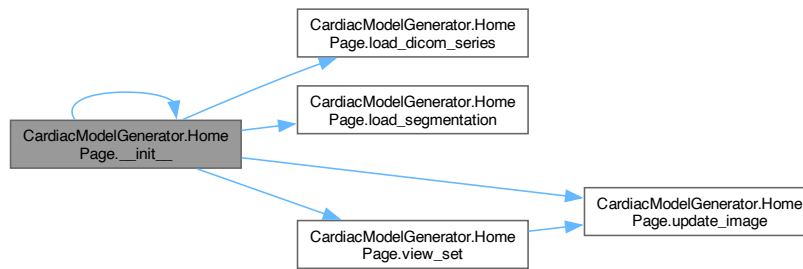

Here is the caller graph for this function:

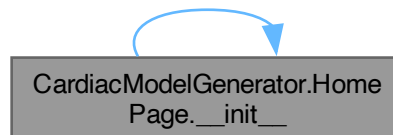

## 5.3.3 Member Function Documentation

### 5.3.3.1 load\_dicom\_series()

```
CardiacModelGenerator.HomePage.load_dicom_series (
    self,
    view_num)
```

**@brief** Loads a DICOM series for a specified view.

**@param** view\_num The view number (1, 2, or 3) to associate with the loaded DICOM series. Example views for heart are 1, 2, and 3.

**@details** Opens a directory dialog to select a folder containing the DICOM series. The series is processed to extract the image volume, coordinates, and individual image objects. These are stored in the parent window's `'dicom_data'` dictionary.

Here is the caller graph for this function:





### 5.3.3.2 load\_segmentation()

```
CardiacModelGenerator.HomePage.load_segmentation (
    self,
    view_num)
```

@brief Loads a segmentation mask for a specified view.

@param view\_num The view number (1, 2, or 3) to associate with the loaded segmentation mask.

@details Opens a file dialog to select a segmentation file. The mask is processed and stored in the parent window's `'segmentation_data'` dictionary.

Here is the caller graph for this function:



### 5.3.3.3 update\_image()


```
CardiacModelGenerator.HomePage.update_image (
    self,
    event)
```

@brief Updates the displayed image based on the slider value.

@param event The wxPython event triggering this action.

@details Extracts the selected slice from the image stack, resizes it to maintain aspect ratio, and updates the image in the `StaticBitmap` widget. The image is optionally rotated for better visualization.

Here is the caller graph for this function:

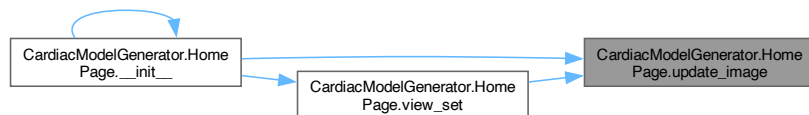

#### 5.3.3.4 view\_set()


```
CardiacModelGenerator.HomePage.view_set (  
    self,  
    view_num)
```

@brief Displays a DICOM image stack with its corresponding segmentation overlay for a specified view.

@param view\_num The view number (1, 2, or 3) to display.

@details Checks if both DICOM images and segmentation masks are loaded for the specified view. If they are, the segmentation is overlaid on the images, and the stack is displayed. The slider is configured to navigate the slices. Displays an error message if the data is missing or invalid.

Here is the call graph for this function:

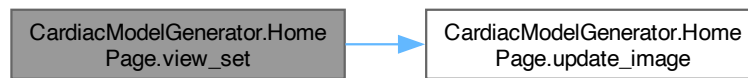

Here is the caller graph for this function:




The documentation for this class was generated from the following file:


- `CardiacModelGenerator.py`

## 5.4 CardiacModelGenerator.PointCloudOptions Class Reference

Inheritance diagram for CardiacModelGenerator.PointCloudOptions:

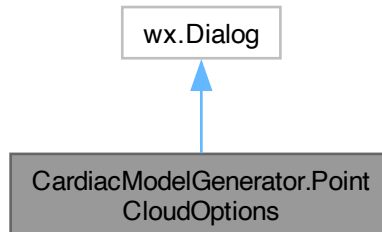

Collaboration diagram for CardiacModelGenerator.PointCloudOptions:

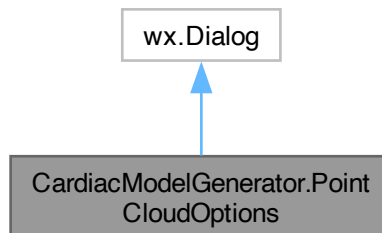

### Public Member Functions

- `__init__` (self, parent, \*args, \*\*kwargs)
- `update_point_size_value` (self, event)
- `update_merging_tolerance_value` (self, event)
- `on_generate_point_cloud` (self, event)

### Public Attributes

- `colormap_combo` = `wx.ComboBox(self, choices=colormap_choices, style=wx.CB_READONLY)`
- `point_size_slider` = `wx.Slider(self, minValue=3, maxValue=20, value=3, style=wx.SL_HORIZONTAL)`
- `point_size_value` = `wx.StaticText(self, label=str(self.point_size_slider.GetValue()))`
- `merging_tolerance_slider`
- `merging_tolerance_value` = `wx.StaticText(self, label="{self.merging_tolerance_slider.GetValue() / 100:.2f}")`
- `whichmask_text` = `wx.TextCtrl(self, value="1")`
- `generate_button` = `wx.Button(self, label="Generate Point Cloud")`

- **on\_generate\_point\_cloud**
- **update\_point\_size\_value**
- **update\_merging\_tolerance\_value**
- **colormap** = None
- **int point\_size** = 3
- **float merging\_tolerance** = 0.0
- **int whichmask** = 1

### 5.4.1 Detailed Description

```
@class PointCloudOptions
@brief Dialog for configuring point cloud generation options.
@details Provides controls for selecting a colormap, adjusting point size, setting merging tolerance,
        and specifying the mask to use for generating the point cloud.

@uml
@startuml
class PointCloudOptions {
    - colormap_combo : wx.ComboBox
    - point_size_slider : wx.Slider
    - point_size_value : wx.StaticText
    - merging_tolerance_slider : wx.Slider
    - merging_tolerance_value : wx.StaticText
    - whichmask_text : wx.TextCtrl
    - generate_button : wx.Button
    - colormap : str
    - point_size : int
    - merging_tolerance : float
    - whichmask : int

    + __init__(parent : wx.Window, *args, **kwargs)
    + update_point_size_value(event : wx.Event)
    + update_merging_tolerance_value(event : wx.Event)
    + on_generate_point_cloud(event : wx.Event)
}

PointCloudOptions -- wx.Dialog : inherits
PointCloudOptions o-- wx.ComboBox : "Colormap selection"
PointCloudOptions o-- wx.Slider : "Adjust point size and merging tolerance"
PointCloudOptions o-- wx.TextCtrl : "Mask selection"
PointCloudOptions o-- wx.Button : "Generate point cloud"
PointCloudOptions --> Colormap : "Uses matplotlib colormap options"
PointCloudOptions --> PointCloud : "Generates point cloud with configured options"

' Notes for context
note top of PointCloudOptions
    PointCloudOptions allows the user to configure parameters for generating a point cloud.
    It includes widgets for selecting colormap, adjusting point size and merging tolerance,
    and specifying the segmentation mask to use.
end note
@enduml
```

### 5.4.2 Constructor & Destructor Documentation

#### 5.4.2.1 \_\_init\_\_()


```
CardiacModelGenerator.PointCloudOptions.__init__ (
    self,
    parent,
    * args,
    ** kwargs)
```

```


@brief Initializes the PointCloudOptions dialog.
@param parent The parent window that contains this dialog.
@param *args Additional positional arguments for the wx.Dialog.
@param **kwargs Additional keyword arguments for the wx.Dialog.
@details Sets up the layout, widgets, and event bindings for configuring point cloud generation options.
         Includes controls for colormap selection, point size adjustment, merging tolerance, and mask selection.

```

Here is the call graph for this function:

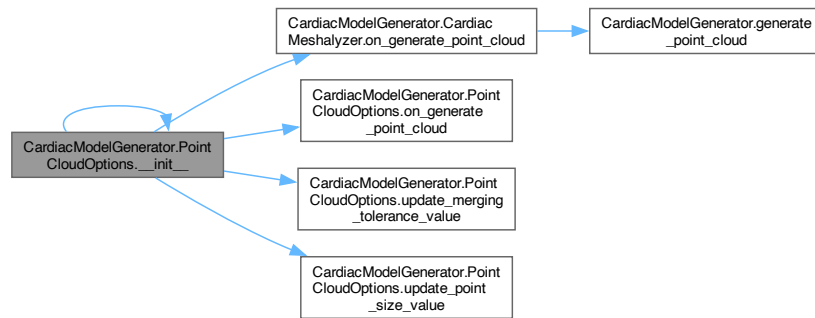

Here is the caller graph for this function:

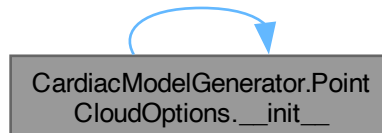

## 5.4.3 Member Function Documentation

### 5.4.3.1 on\_generate\_point\_cloud()

```

CardiacModelGenerator.PointCloudOptions.on_generate_point_cloud (
    self,
    event)


```

```

@brief Retrieves user inputs and closes the dialog.
@param event The wxPython button event triggering this action.
@details Captures user-selected values from the colormap dropdown, point size slider, merging tolerance slider
         and mask text input. Validates the mask input as an integer and closes the dialog with a success status.
@throws ValueError If the mask input is not a valid integer.

```

Here is the caller graph for this function:

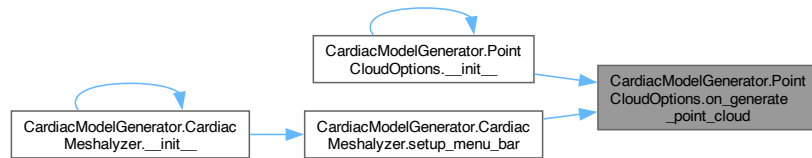

#### 5.4.3.2 update\_merging\_tolerance\_value()


```
CardiacModelGenerator.PointCloudOptions.update_merging_tolerance_value (
    self,
    event)
```

@brief Updates the displayed value for the merging tolerance slider.

@param event The wxPython slider event triggering this action.

@details Converts the slider value to a floating-point representation (0.00 to 5.00) and updates the 'merging\_tolerance\_value' label for user feedback.

Here is the caller graph for this function:

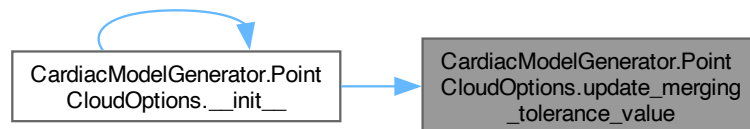

#### 5.4.3.3 update\_point\_size\_value()


```
CardiacModelGenerator.PointCloudOptions.update_point_size_value (
    self,
    event)
```

@brief Updates the displayed value for the point size slider.

@param event The wxPython slider event triggering this action.

@details Reflects the current slider value in the 'point\_size\_value' label to provide immediate feedback to the user.

Here is the caller graph for this function:

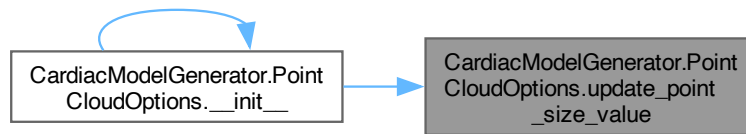

## 5.4.4 Member Data Documentation

### 5.4.4.1 merging\_tolerance\_slider

`CardiacModelGenerator.PointCloudOptions.merging_tolerance_slider`

#### Initial value:


```
= wx.Slider(  
    self, minValue=0, maxValue=500, value=0, style=wx.SL_HORIZONTAL  
)
```

The documentation for this class was generated from the following file:


- `CardiacModelGenerator.py`

## 5.5 CardiacModelGenerator.StartPage Class Reference

Inheritance diagram for `CardiacModelGenerator.StartPage`:

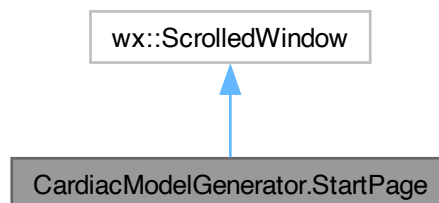

Collaboration diagram for CardiacModelGenerator.StartPage:

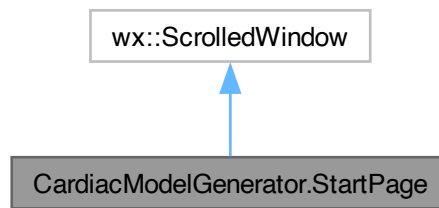

### Public Member Functions

- `__init__` (self, parent)
- `load_image` (self)
- `on_resize` (self, event)
- `open_home_page` (self, event)
- `close_program` (self, event)

### Public Attributes

- `sizer` = `wx.BoxSizer(wx.VERTICAL)`
- `image_path` = `IMAGE_PATH`
- `bitmap` = `wx.StaticBitmap(self)`
- `close_program`
- `open_home_page`
- `on_resize`

## 5.5.1 Detailed Description

```

@class StartPage
@brief Introductory page for the Cardiac Meshalyzer application.
@details This class represents the starting page of the GUI, which includes a title, an introductory image,
        a warning message, and navigation buttons for proceeding or exiting the application.

@uml
@startuml
class StartPage {
    - sizer : wx.BoxSizer
    - image_path : str
    - bitmap : wx.StaticBitmap

    + __init__(parent : wx.Window)
    + load_image()
    + on_resize(event : wx.Event)
    + open_home_page(event : wx.Event)
    + close_program(event : wx.Event)
}

' Associations to other elements in the GUI
StartPage *-- wx.ScrolledWindow : inherits
StartPage o-- wx.BoxSizer : "Main vertical sizer"
  
```



```

StartPage o-- wx.StaticBitmap : "Image placeholder"
StartPage --> wx.Button : "Handles Close and Continue buttons"

' Notes for additional context
note top of StartPage
    StartPage serves as the introductory page for the Cardiac Meshalyzer GUI.
    It includes a title, image placeholder, warning text, and navigation buttons.
    The buttons allow users to navigate to the Home Page or close the application.
end note

@enduml

```

## 5.5.2 Constructor & Destructor Documentation

### 5.5.2.1 \_\_init\_\_()


```

CardiacModelGenerator.StartPage.__init__ (
    self,
    parent)


@brief Initializes the StartPage class and its UI components.
@param parent The parent wx.Window to which this StartPage belongs.
@details This method sets up the introductory page of the Cardiac Meshalyzer application, including:
    - A black background with a title.
    - A placeholder for an introductory image.
    - A warning message about using the application.
    - Two navigation buttons ("Close" and "Continue").
    Enables scrolling and binds events for resize, navigation, and program termination.

```

Here is the call graph for this function:

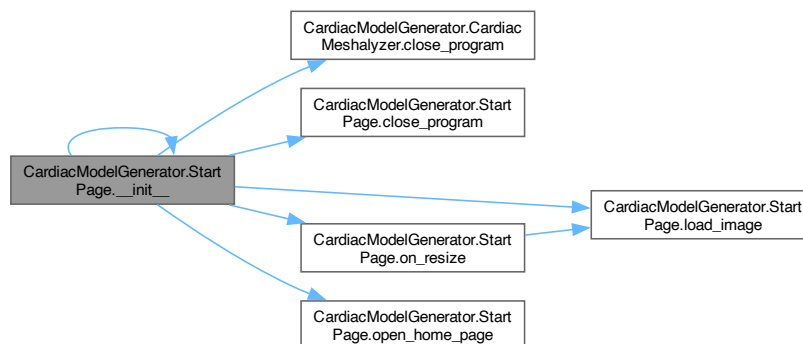

Here is the caller graph for this function:

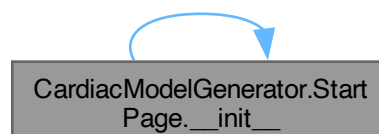


## 5.5.3 Member Function Documentation

### 5.5.3.1 close\_program()

```
CardiacModelGenerator.StartPage.close_program (
    self,
    event)
```

@brief Closes the application.  
 @param event The wxPython event triggering this action.  
 @details Calls the 'Close' method of the top-level window to terminate the application.

Here is the caller graph for this function:

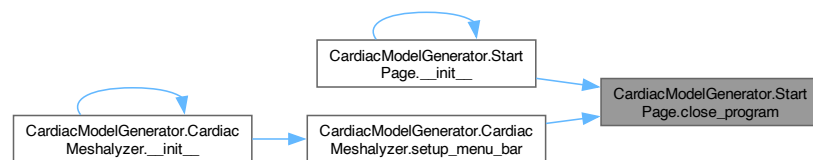


### 5.5.3.2 load\_image()

```
CardiacModelGenerator.StartPage.load_image (
    self)
```

@brief Loads and displays the introductory image.  
 @details Attempts to load the image from the specified path and scales it to fit the current window size.  
 If the image cannot be loaded, a placeholder message is displayed.  
 @throws Exception If the image file cannot be read or scaled.

Here is the caller graph for this function:

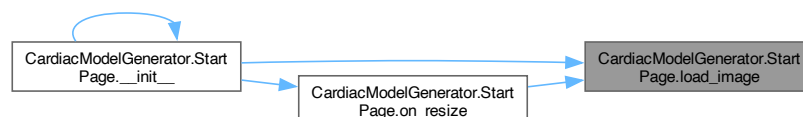

### 5.5.3.3 on\_resize()

```
CardiacModelGenerator.StartPage.on_resize (  
    self,  
    event)
```

@brief Handles window resize events.


@param event The wxPython resize event triggering this action.

@details Reloads and rescales the introductory image to fit the updated window size.

Here is the call graph for this function:



Here is the caller graph for this function:

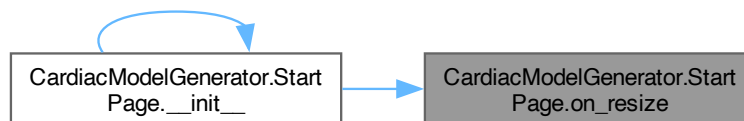

### 5.5.3.4 open\_home\_page()


```
CardiacModelGenerator.StartPage.open_home_page (  
    self,  
    event)
```

@brief Navigates to the Home Page of the application.

@param event The wxPython event triggering this action.

@details Calls the 'show\_page' method of the top-level window to display the "Home Page".

Here is the caller graph for this function:

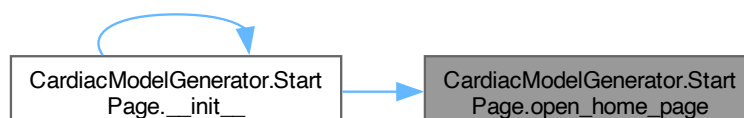

The documentation for this class was generated from the following file:

- `CardiacModelGenerator.py`

# Index

- `__init__`
    - `CardiacModelGenerator.CardiacMeshalyzer`, [13](#)
    - `CardiacModelGenerator.CleanTetraMeshOptions`, [24](#)
    - `CardiacModelGenerator.HomePage`, [27](#)
    - `CardiacModelGenerator.PointCloudOptions`, [32](#)
    - `CardiacModelGenerator.StartPage`, [37](#)
  - `_init_`, [7](#)
- `add_page`
  - `CardiacModelGenerator.CardiacMeshalyzer`, [14](#)
- `CardiacModelGenerator`, [7](#)
  - `clean_tetra_mesh`, [8](#)
  - `generate_point_cloud`, [8](#)
  - `generate_tetra_mesh`, [9](#)
  - `get_cell_quality`, [9](#)
  - `save_vista`, [9](#)
- `CardiacModelGenerator.CardiacMeshalyzer`, [11](#)
  - `__init__`, [13](#)
  - `add_page`, [14](#)
  - `clear_all_files`, [15](#)
  - `close_program`, [15](#)
  - `get_masks`, [15](#)
  - `getMaskOverlay`, [16](#)
  - `ImportDicomSeries`, [16](#)
  - `on_clean_tetra_mesh`, [16](#)
  - `on_extract_mesh_quality`, [17](#)
  - `on_generate_point_cloud`, [17](#)
  - `on_generate_tetra_mesh`, [18](#)
  - `open_point_cloud_options`, [19](#)
  - `save_point_cloud`, [19](#)
  - `save_tetra_cloud`, [19](#)
  - `setup_menu_bar`, [20](#)
  - `show_page`, [21](#)
- `CardiacModelGenerator.CleanTetraMeshOptions`, [22](#)
  - `__init__`, [24](#)
  - `on_cancel`, [25](#)
  - `on_ok`, [25](#)
- `CardiacModelGenerator.HomePage`, [26](#)
  - `__init__`, [27](#)
  - `load_dicom_series`, [28](#)
  - `load_segmentation`, [28](#)
  - `update_image`, [29](#)
  - `view_set`, [29](#)
- `CardiacModelGenerator.PointCloudOptions`, [31](#)
  - `__init__`, [32](#)
  - `merging_tolerance_slider`, [35](#)
  - `on_generate_point_cloud`, [33](#)
  - `update_merging_tolerance_value`, [34](#)
  - `update_point_size_value`, [34](#)
- `CardiacModelGenerator.StartPage`, [35](#)
  - `__init__`, [37](#)
  - `close_program`, [38](#)
  - `load_image`, [38](#)
  - `on_resize`, [38](#)
  - `open_home_page`, [39](#)
- `clean_tetra_mesh`
  - `CardiacModelGenerator`, [8](#)
- `clear_all_files`
  - `CardiacModelGenerator.CardiacMeshalyzer`, [15](#)
- `close_program`
  - `CardiacModelGenerator.CardiacMeshalyzer`, [15](#)
  - `CardiacModelGenerator.StartPage`, [38](#)
- `generate_point_cloud`
  - `CardiacModelGenerator`, [8](#)
- `generate_tetra_mesh`
  - `CardiacModelGenerator`, [9](#)
- `get_cell_quality`
  - `CardiacModelGenerator`, [9](#)
- `get_masks`
  - `CardiacModelGenerator.CardiacMeshalyzer`, [15](#)
- `getMaskOverlay`
  - `CardiacModelGenerator.CardiacMeshalyzer`, [16](#)
- `ImportDicomSeries`
  - `CardiacModelGenerator.CardiacMeshalyzer`, [16](#)
- `load_dicom_series`
  - `CardiacModelGenerator.HomePage`, [28](#)
- `load_image`
  - `CardiacModelGenerator.StartPage`, [38](#)
- `load_segmentation`
  - `CardiacModelGenerator.HomePage`, [28](#)
- `merging_tolerance_slider`
  - `CardiacModelGenerator.PointCloudOptions`, [35](#)
- `on_cancel`
  - `CardiacModelGenerator.CleanTetraMeshOptions`, [25](#)
- `on_clean_tetra_mesh`
  - `CardiacModelGenerator.CardiacMeshalyzer`, [16](#)
- `on_extract_mesh_quality`
  - `CardiacModelGenerator.CardiacMeshalyzer`, [17](#)
- `on_generate_point_cloud`
  - `CardiacModelGenerator.CardiacMeshalyzer`, [17](#)
  - `CardiacModelGenerator.PointCloudOptions`, [33](#)
- `on_generate_tetra_mesh`
  - `CardiacModelGenerator.CardiacMeshalyzer`, [18](#)

on\_ok  
    CardiacModelGenerator.CleanTetraMeshOptions,  
        [25](#)

on\_resize  
    CardiacModelGenerator.StartPage, [38](#)

open\_home\_page  
    CardiacModelGenerator.StartPage, [39](#)

open\_point\_cloud\_options  
    CardiacModelGenerator.CardiacMeshalyzer, [19](#)

save\_point\_cloud  
    CardiacModelGenerator.CardiacMeshalyzer, [19](#)

save\_tetra\_cloud  
    CardiacModelGenerator.CardiacMeshalyzer, [19](#)

save\_vista  
    CardiacModelGenerator, [9](#)

setup\_menu\_bar  
    CardiacModelGenerator.CardiacMeshalyzer, [20](#)

show\_page  
    CardiacModelGenerator.CardiacMeshalyzer, [21](#)

update\_image  
    CardiacModelGenerator.HomePage, [29](#)

update\_merging\_tolerance\_value  
    CardiacModelGenerator.PointCloudOptions, [34](#)

update\_point\_size\_value  
    CardiacModelGenerator.PointCloudOptions, [34](#)

view\_set  
    CardiacModelGenerator.HomePage, [29](#)